



PHOTON
SURGICAL
SYSTEMS

THE ECOMAX X-RAY FILM PROCESSOR



**Ecological and Economical
Film Processor**

THE ECOMAX X-RAY FILM PROCESSOR


FEATURES

- Compact size, ideal for small darkrooms
- Aesthetically pleasing and technically innovative
- Designed to offer optimal image quality
- Clever tank design reduces chemical oxidation
- Ingeniously simple to use and maintain

SPECIFICATIONS

- Processing Time 100 seconds (factory set)
- Film Sizes 10x10cm to 35x43cm
- Capacity 74 films per hour (24x30cm)
- Film Types Any film suitable for 100 second processing
- Chemistry Any chemistry suitable for rapid processing
- Water consumption 1.5 litres per minute (during processing)
- Dimensions Length: 83.8cm Width: 66cm Height: 42cm
- Weight 26.5 kg (empty)

Specification subject to change without prior notice (April 2010)

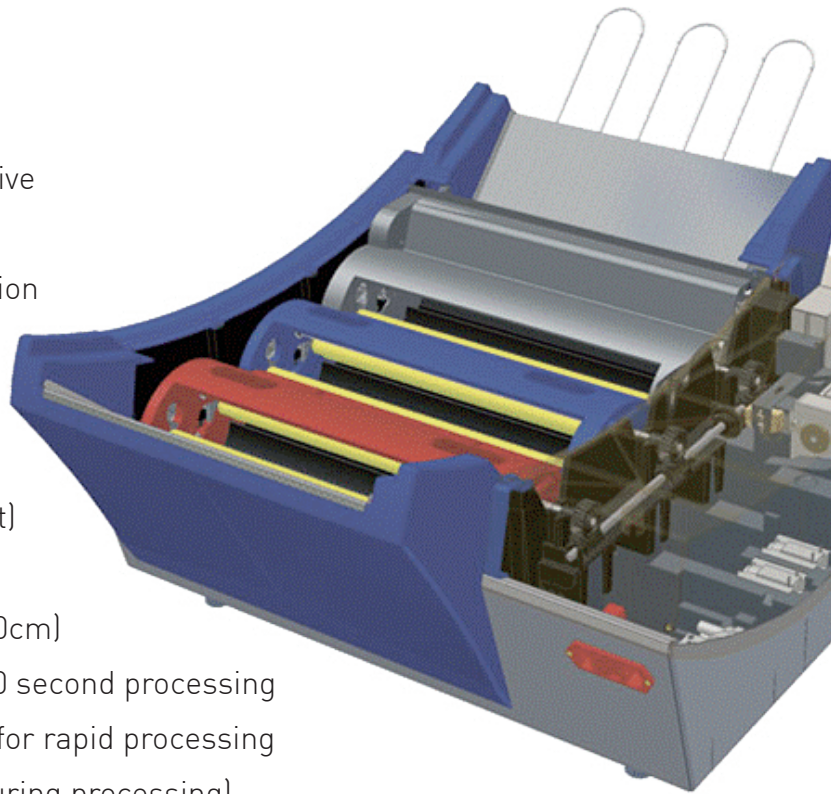


BENEFITS

- All components are easily accessible reducing down time and maintenance costs
- Innovative roller rack ensures a more gentle film path
- Combination of innovative roller rack design and newly shaped chemical tanks result in optimal image quality

ECOLOGICAL AND ECONOMICAL

- Smaller chemical tanks reduce chemistry oxidation of large quantities of chemistry saving money and making it more environmentally friendly
- Anti-crystallisation, anti-oxidation and permanent circulation of the chemistry functions included
- Automatic replenishment adjustment
- Rests in stand-by mode when not in use

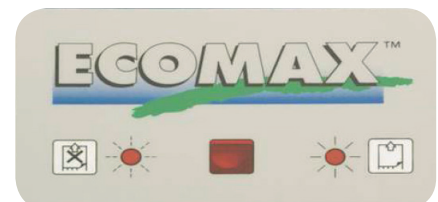

Acoustic sound and light display when ready to process film.



Telephone: 01452 347101

E-mail: info@photonsurgicalsystems.co.uk

Aquarius Centre, Edison Close, Waterwells Business Park, Quedgeley, Gloucestershire, GL2 2FN